

DusQ 1 CPG 1000

<http://www.lumiprobe.com/p/dusq-1-cpg-1000>


This support with a pore size of 1000 Å is intended for the automated synthesis of oligonucleotides of up to 120 bases in length modified with DusQ1 dark quencher at the 3' end.

Non-fluorescent DusQ1 quencher exhibits the strongest absorption within the range of 480 to 580 nm; its absorption maximum is at 534 nm. It is suitable for quenching (a combination of static and dynamic quenching) of many fluorophores, including Biosearch Blue™, Marina Blue™, Edans, Bothell Blue, FAM™, JOE™, VIC™, R6G, HEX™, TET™ and Yakima Yellow™. It can be used for the synthesis of hybridization probes such as TaqMan, Molecular Beacon, Scorpion.


Usage

Coupling: Standard conditions identical to normal nucleobases.

Deprotection: 2 hours at room temperature using concentrated ammonia or 10 min at 65 °C using AMA mixture, concentrated aqueous ammonia/40% methylamine (1:1). Deprotection conditions depend on oligonucleotide composition and nucleobase protecting groups, as well as additional modifications, if present.

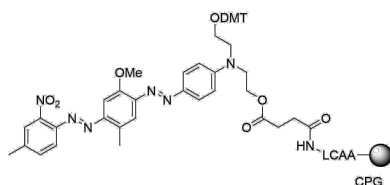

Structure of DusQ 1 CPG 1000

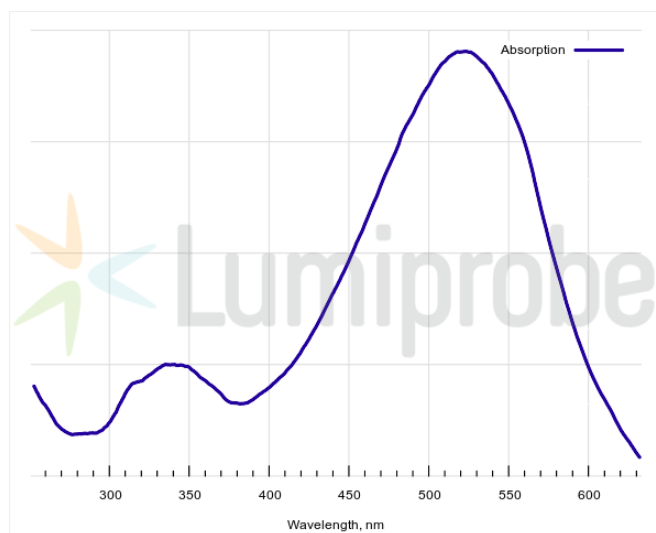

Absorption spectrum of DusQ 1

General properties

Appearance:	red beads
Quality control:	NMR ¹ H and HPLC-MS (95%) of bound reagent, loading measurement, functional testing in oligo synthesis.
Storage conditions:	Storage: 24 months after receipt at -20°C in the dark. Transportation: at room temperature for up to 3 weeks. Avoid prolonged exposure to light. Desiccate.

Spectral properties

Excitation/absorption maximum, nm:	522
ε, L·mol ⁻¹ ·cm ⁻¹ :	27300
CF ₂₆₀ :	0.17
CF ₂₈₀ :	0.10

Oligo synthesis details

Pore size, Å:	1000
---------------	------

Typical loading, $\mu\text{mol/g}$:

30–50